


タンパク質の分離(2)

SunShell C4-30 2.6 μ m, 100 x 2.1 mm i.d.

Protein Mix using SunShell C4-30

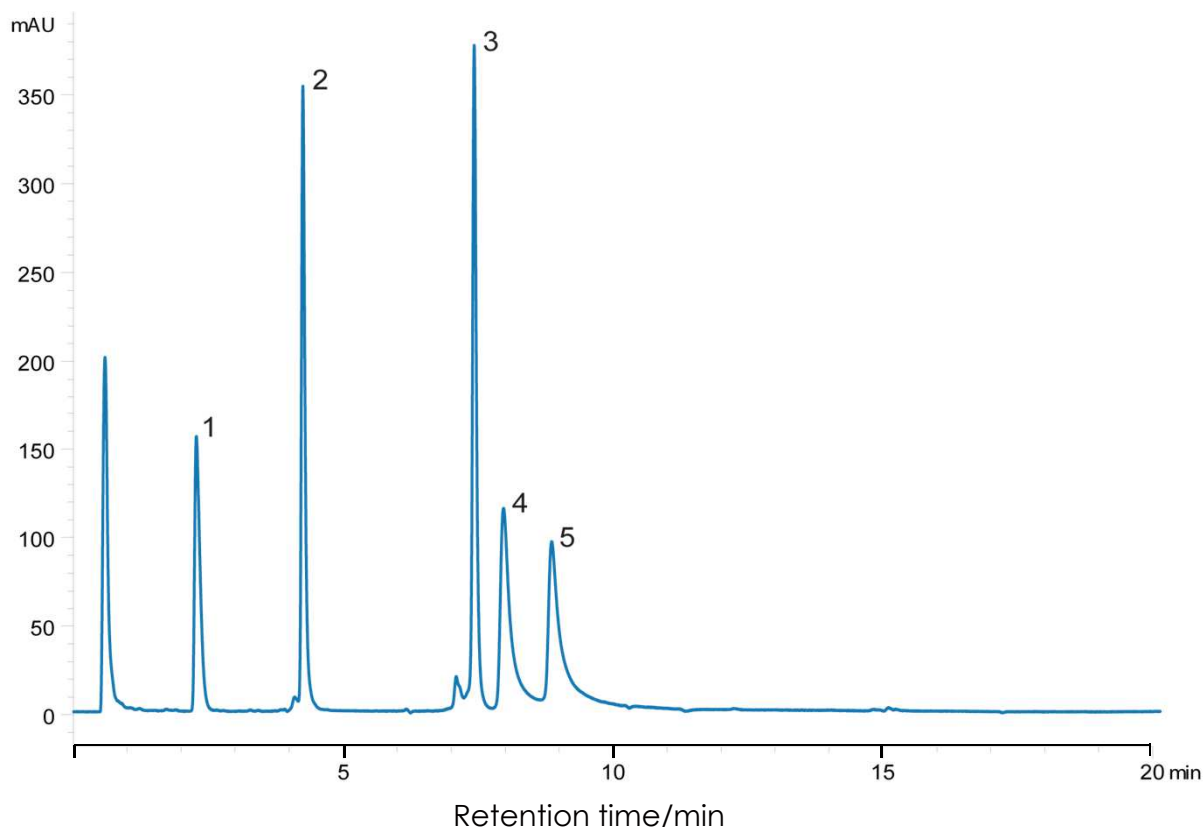

Column: SunShell C4-30 2.6 μ m (30 nm), 100 x 2.1 mm

Mobile phase: A) 0.1% TFA in water

B) 0.1% TFA in acetonitrile

Gradient program:	Time(min)	0	20
	%B	22	70.5

Flow rate: 0.50 mL/min

Temperature: 40 °C

Detection: UV@214 nm

Injection volume: 10 μ L

Concentration: 0.01 μ g/ μ L each protein

Sample:

1. Angiotensin I
2. Ribonuclease A
3. Lysozyme
4. Transferrin
5. Bovine Serum Albumin